


1 - Open Oculus App, make sure the device is connected and set up properly. Enable “Oculus link” if you are using a Quest/Quest2;

2 - Open “Anaconda Prompt”, type:

- `conda activate psychxr`
- `spyder`

3 - Run “fold_parallax_v23.py” script, a dialogue will show up on the task bar (flashing orange, click on it):

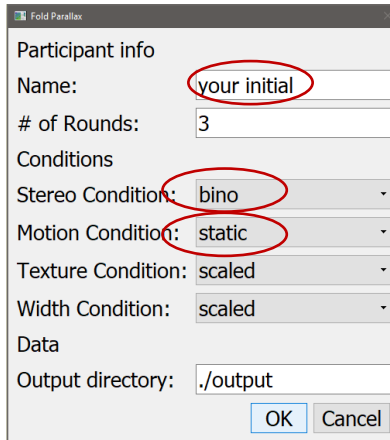

Participant info

Name: your initial

of Rounds: 3

Conditions

Stereo Condition: bino

Motion Condition: static

Texture Condition: scaled

Width Condition: scaled

Data

Output directory: ./output

OK Cancel

4 - Put 2 letters in “Name” field;

5 - Vary “Stereo Condition” and “Static Condition” (4 conditions in total).